

3052-GEM STEREO GEM STAND MICROSCOPE

**The Portable
Microscope
For the
Gemologist
On-the-Go**

COMPACT
DESIGN
FOR EASY
TRANSPORT

SIMPLE
TO USE

GREAT
VALUE
AT A
GREAT
PRICE!



Delivering quality through vision

3052-GEM STEREO GEM STAND MICROSCOPE

ACCU-SCOPE's 3052-GEM Stereo Gem Stand Microscope is compact, easily transportable and is ideal for the gemologist on the go!

The 3052-GEM is a simple, cost-effective solution for the discerning gemologist looking for a great value at an affordable price.

Easily observe loose stones with the included wire stone holder, or view rings with optional the ring holder (sold separately).



Accessories shown left:

- ① Inclusion pointer (optional)
- ② Ring holder (optional)
- ③ Wire stone holder (included)


SPECIFICATIONS

Optical System	1x/3x fixed magnifications for viewing stones at 10x and 30x magnification with the included WF 10x/20mm eyepieces Magnification can be increased to 15x and 45x with the optional 15x eyepieces (sold separately)
Illumination	Built-in, rechargeable LED illumination provides darkfield illumination and reflected illumination from above
Portability/Stand	Powered by a light weight 110v or optional 220v (sold separately) external power supply and rechargeable AA batteries, allows the user to work anywhere
Accessories	Steel wire stone holder allows for multi-angled viewing of various sized gems is included; optional ring holder or inclusion pointer are also available and sold separately
Private Label	This microscope may be private labeled upon request
Warranty	5 year limited warranty for parts and labor, and electronic components

Design, features and specifications are subject to change without notice.

The included wire Stone holder swings out for easier changing of gemstones and rotates for viewing at any angle.

An optional inclusion pointer is also available.

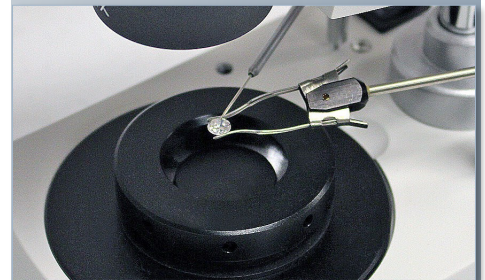

The optional ring holder also pivots & rotates 360° for viewing at any angle.



ISO 9001

Design and production adheres to ISO 9001 international quality standard.

ISO 14001

Design and production meets the Requirements of international standard ISO 14001 for environmental management.



73 Mall Drive • Commack, NY 11725
631-864-1000 (P) • 631-543-8900 (F)
email: info@accu-scope.com • www.accu-scope.com